



M13 bacteriophage (g3p) Biotinylated. Mouse Monoclonal Antibody

Coat protein A; g3p; Minor coat protein

BACKGROUND

The display of repertoires of antibody fragments on the surface of filamentous phage offers a new way to produce immunoreagents with defined specificities.


Phage derived antibody fragments offer a number of advantages over mouse monoclonal antibodies, such as better clearance from the blood, the possibility to select from human combinatorial libraries and the relative ease by which such fragments can be manipulated. The phage display technique thus facilitates the selection of antibody fragments of therapeutic value or research interest.

Antibodies to M13 filamentous phage coat proteins are instrumental in the selection and detection of phages expressing specific antibody fragments or peptide sequences at their surface.

IMMUNOGEN

Hybridoma produced by the fusion of splenocytes from mice immunized with isolated M13 phage coat proteins and mouse myeloma cells.

Data represents absorbancy readings for A10B phage on rabbit IgG (A10B/IgG), A10B phage on BSA (A10B/BSA), streptavidin on rabbit IgG (SA/IgG) and streptavidin on BSA (SA/BSA) for each dilution of biotinylated anti-M13 monoclonal antibody.

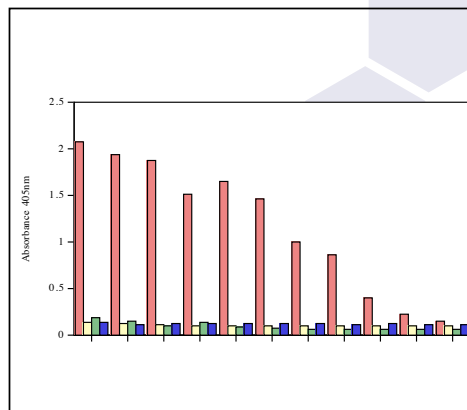

ORDERING INFORMATION

CATALOG NUMBER
Z110M

SIZE
100 µg
FORM
Biotin

HOST/CLONE
Mouse Clone E1

FORMULATION
Provided as solution in phosphate buffered saline with 0.08% sodium azide

CONCENTRATION
See vial for concentration

ISOTYPE
IgG2a

APPLICATIONS
ELISA, Western Blot, Flow Cytometry, IHC

SPECIES REACTIVITY
N/A

ACCESSION NUMBER
M13 Phage P69168

POSITIVE CONTROL/TISSUE EXPRESSION

COMMENTS

Antibody specific for the g3p protein. Antibody can be used for immunohistochemistry, Western blot (1-5 µg/ml), Flow cytometry (1 µg/10⁶ cells) and ELISA. Optimal concentration should be evaluated by serial dilutions.

PURIFICATION

Protein A/G Chromatography

SHIP CONDITIONS

Ship at ambient temperature, freeze upon arrival

STORAGE CUSTOMER

Product should be stored at -20°C. Aliquot to avoid freeze/thaw cycles

STABILITY

Products are stable for one year from purchase when stored properly

REFERENCES

- 1- Van Wezenbeek P.M., Schoenmakers J.G.; Nucleotide sequence of the genes III, VI and I of bacteriophage M13; Nucleic Acids Res. 6:2799-2818(1979).
- 2- Van Wezenbeek P.M., Hulsebos T.J., Schoenmakers J.G.; Nucleotide sequence of the filamentous bacteriophage M13 genome: comparison with phage fd; Gene 11:129-148(1980).
- 3- Cleary J.M., Ray D.S.; Deletion analysis of the cloned replication origin region from bacteriophage M13; J. Virol. 40:197-203(1981).
- 4- Hines J.C., Ray D.S.; Construction and characterization of new coliphage M13 cloning vectors; Gene 11:207-218 (1980).
- 5- Yanisch-Perron C., Vieira J., Messing J.; Improved M13 phage cloning vectors and host strains: nucleotide sequences of the M13mp18 and pUC19 vectors; Gene 33:103-119(1985).
- 6- Messing J.; New M13 Vectors for Cloning; Meth. Enzymol. 101:20 78(1983).
- 7- Sanger F., Nicklen S., Coulson A.R.; DNA sequencing with chain-terminating inhibitors; Proc. Natl. Acad. Sci. U.S.A. 74:5463-5468(1977).
- 8- Zoller M.J., Smith M.; Oligonucleotide-directed mutagenesis of DNA fragments cloned into M13 vectors; Meth. Enzymol. 100:468-500(1983).
- 9- Hu N.T., Messing J.; The making of strand-specific M13 probes; Gene 17:271-277(1982).
- 10- Heidecker G., Messing J., Gronenborn B.; A versatile primer for DNA sequencing in the M13mp2 cloning system; Gene 10:69-73(1980).
- 11- Ebricht R., Dong Q., Messing J.; Corrected nucleotide sequence of M13mp18 gene III; Gene 114:81-83(1992).
- 12- Hong G.F.; A method for sequencing single-stranded cloned DNA in both directions; Biosci. Rep. 1:243-252(1981).
- 13- Messing J.; Multipurpose cloning system based on the single-stranded DNA bacteriophage M13; Recombinant DNA Technical Bulletin NIH 2:43-48(1979).

PRODUCT SPECIFIC REFERENCES